

# Rocket<sup>®</sup> IPC<sup>™</sup> Drainage Kit

## INSTRUCTIONS FOR USE - PREPARATION & DRAINAGE



1

Open bottle pack, place bottle on edge of wrap, ensure tubing stays on blue drape



Ready



Do NOT use

Check bottle indicator, if raised (above), use a new bottle



3

Put on the sterile gloves, handling the cuffs only



4

a. Peel open the valve cover pack.  
b. Tear open 3 alcohol wipes.  
c. Remove clear dressing from pack. Place items into centre of blue wrap



5

Remove tip cover and place exposed tip on blue wrap




6

Remove valve cap



7

Clean with new alcohol swab

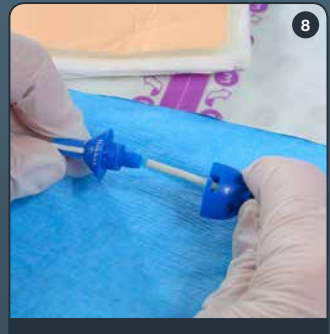

8

Gently insert access tip into valve, twist to lock

# Rocket<sup>®</sup> IPC<sup>™</sup> Drainage Kit


## INSTRUCTIONS FOR USE - PREPARATION & DRAINAGE



Slide white clamp to open position



Push button to drain

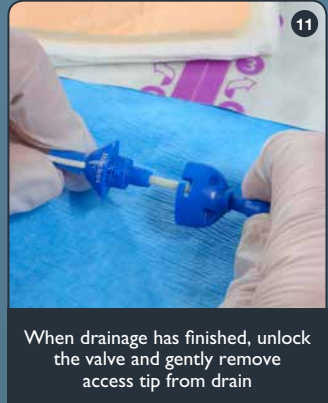


When drainage has finished, unlock the valve and gently remove access tip from drain



Using an alcohol swab, gently clean around the outside of the valve



Replace new valve cap

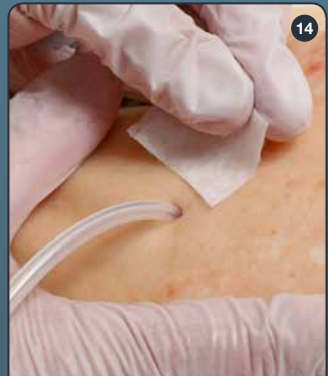

Clean wound site carefully with an alcohol swab



Coil the drain over the foam dressing



Cover coiled drain with gauze squares and place clear dressing over them



Dressing applied