

Rocket® Spinal Manometer Sets

Diagnostic Medicine

with **NRFit™** connectors

ROCKET® SPINAL MANOMETERS with NRFit™ ensure connection to ISO80369-6 compliant spinal needles and syringes for the safe removal and measurement of CSF.

STACKABLE ASSEMBLY: NRFit™ Male/female fittings of each tube permit additional sections to be stacked to allow pressures $>40\text{cmH}_2\text{O}$ to be measured when required.

CLEAR, STURDY PLASTIC CONSTRUCTION for reliable, secure assembly and ease of use.


SUPPLIED COMPLETE with 3 way stopcock for attachment to all NRFit™ ISO80369-6 compliant spinal needles and syringes.

READY FOR IMMEDIATE USE for when time is precious.

Rocket® Spinal Manometers have been providing safe, reliable, cost effective CSF pressure measurement for over 30 years.

The unique 'stackable' design, original to Rocket, has become a standard neurological diagnostic tool and can be found in Emergency Departments, Neurological, Medical Wards and specialist Neurosurgery centres throughout the world.

This original Rocket® design is now available with ISO80369-6 NRFit™ compliant connectors.

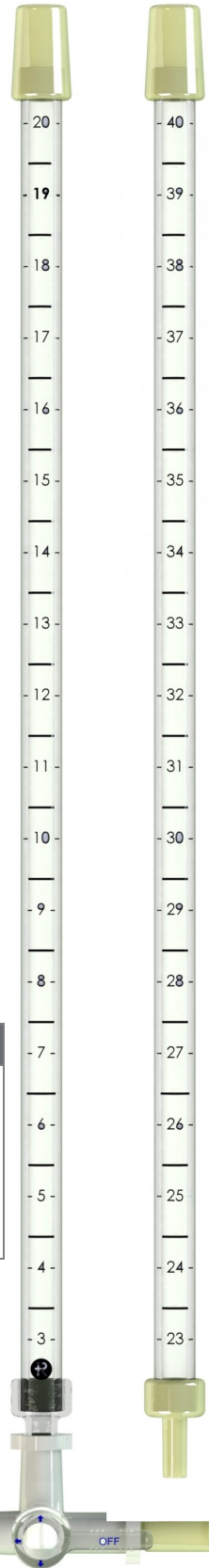

Description	Order Code	NHSSC Code
Rocket® Spinal Manometer Set with NRFit™ connectors Sterile, for single use, packed in cartons 10 units	R55990-ISO-6	eDirect GBR000

NRFit™ is a trademark of GEDSA and used with permission.

Rocket Medical Plc

Sedling Road, Washington, Tyne & Wear, NE38 9BZ

T: +44 (0) 191 419 6988 E: customerservices@rocketmedical.com W: www.rocketmedical.com

Rev: 06
Date: 15/05/2020
Copyright© 2017-20