


Rocket KCH™ Fetal Bladder Drain

The Rocket® KCH™ Catheter is intended for use in fetal urinary tract decompression following diagnosis of fetal obstructive uropathy in foetuses of 18-38 weeks gestation.



- **TWIN PIGTAIL** catheter with 6 holes for reliable drainage and metal tip markers for easier post insertion position monitoring.
- **MRI COMPATIBLE:** The catheter and marker rings are compatible with MR Imaging for simplified post insertion monitoring.
- **PTFE GUIDEWIRE** supplied to ease loading of pigtail catcher into introducer and minimises the risk of kinks
- **REUSABLE INTRODUCER SET** for secure, insertion and reliable placement.

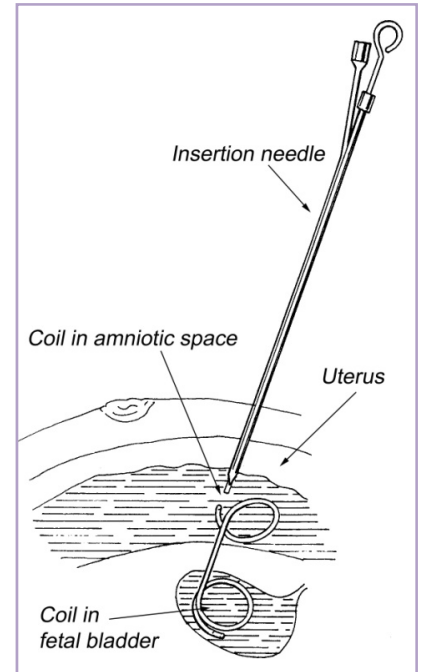

Description	Rocket Product Code	CH2 Order Code
Rocket® KCH™ Fetal Bladder Drain Set with 0.035" x 80cm PTFE guidewire, supplied sterile, pouch packed, each	R57405	1969187
Rocket® KCH™ Reusable Introducer Set: trocar, cannula, long & short pusher rods	R29100	



CLIFFORDHALLAMHEALTHCARE

1300 720 274 www.ch2.net.au



Rocket Medical Pty. Ltd. Suite 209. | Katherine Street, Chatswood, NSW. 2067
T: +61 431 146448 **E:** aus@rocketmedical.com **www:** rocketmedical.com

Company Nos: ABN 65150984937 ACN 50984937