




- ◆ **LX<sup>™</sup> LONGITUDINAL ECHO BEVEL DESIGN** with in-line high resolution echogenic surface preparation to give smoother, easier tissue penetration with pinpoint ultrasound definition of the needle tip.
- ◆ **HUB MARKINGS** ensure correct orientation of the inner stylet following removal and replacement at any time during the procedure.
- ◆ **SECURE LUER LOCK FITTINGS** ensure needle, extension set and syringes remain correctly and securely connected throughout.
- ◆ **INTEGRAL TUBING SET** – purpose designed extension tube set is included in every pack for simple, safe attachment to the needle without the need to find costly off-the-shelf items which cannot offer guaranteed compatibility and security.

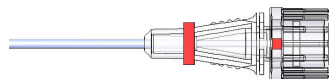

The Rocket<sup>®</sup> LX<sup>™</sup> CVS Needle Set has been developed with the assistance of a leading UK Fetal Medicine Unit and is the result of a 4 year program of continuous product development and design improvement.

**Indications:** For the aspiration of fluid from the amniotic sac and transabdominal chorionic villus sampling for the detection of chromosomal abnormalities in the fetus.



Description	Order Code
<b>Rocket<sup>®</sup> LX<sup>™</sup> Chorionic Villus Sampling Set</b> 17G x 17cm linear echo CVS needle, 10cm luer lock extension tube set. Sterile for single use, packed in cartons of 10 sets	<b>R57034</b>