

Customer Information Bulletin


Rocket Medical® Ambulatory bag

Dear Customer,

Due to ongoing supply chain issues, you may receive the Rocket Medical® R54564 Ambulatory Bag without an activated carbon filter.

The activated carbon filter is there to absorb possible odours from the pleural fluid or empyema and does not compromise the performance of the device with regards to air leak resolution.

You will see that underneath the activated carbon filter is the actual oleophobic filter material which allows air to pass out of the device and the Ambulatory bag to perform as usual.

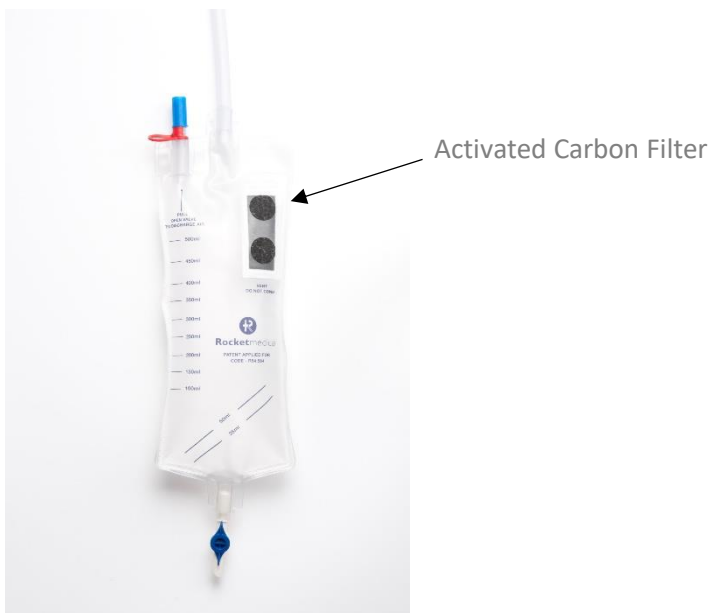

If you have any questions, please contact the Customer Services Team for assistance:

T: +44 (0)191 419 6988 E: customerservices@rocketmedical.com

Customer Information Bulletins for all of our products can be found on our website: www.rocketmedical.com

TEMPLATE153 Rev1 22/FEB/2022