

# Rocket Water Trap Set


## for use with **Craft™** Suction Pumps

The Rocket Water Trap Set for use with **Craft™** Suction Pumps offers a simple, cost effective solution to contamination risks and reuse problems.

- **COMPLETE SET:** Provides the complete bottle and tubing set in single, sterile pack for one time patient use.
- **GLASS TRAP BOTTLE:** To maintain clear view in the event of potential pump contamination with aspirate. Designed to withstand >450mmHg vacuum.
- **REMOVES CONTAMINATION RISK** from irregular or uncontrolled decontamination of reusable bottle, bung and angle tube set.
- **PATIENT TUBE SET:** 2m patient tube set with integral luer lock connector to ensure secure, reliable connection to the aspiration needle set.
- **SUITABLE FOR USE:** Compatible with all versions of **Craft™ Suction and Duo Vac.**  
Also compatible with the *genX and IVF On-line Pioneer Pro Pumps.*



- Disconnect silicone tubing set ❶ from vacuum port
- Discard bung, glass bottle set and tubing ❷, ❸ & ❹
- Replace with water trap set and reconnect tubing at ❶

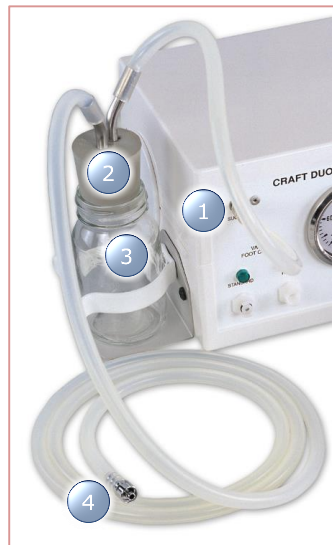

### Ordering Information

**Rocket Water Trap Set:**  
suitable for all **Craft™** Suction Pumps Sterile, for single patient use.

10 per carton

**R57685**

