


# Rocket

## HSG, HyCoSy & Sonohystogram Catheters

The Rocket HSG, HyCoSy & Sonohystogram Catheter range provides a reliable, comfortable and cost effective method of performing uterine canalisation for hysterosalpingography and hysterosalpingo-contrast sonography.

Rocket HSG, HyCoSy & Sonohystogram Catheters have these well-established and trusted features:

- **LATEX FREE** balloon minimises risk of patient tissue reactions.
- **REMOVABLE STIFFENING SHEATH** to aid insertion in the heavily retroverted or antverted cervix.
- **CONTRAST PATH TAP** ensures retention of contrast when further volumes of media are required.
- **SMOOTH TIP & LOW PROFILE BALLOON** ensures low trauma insertion and removal – ideal for the nulliparous cervix.
- **INTEGRAL WIRE STYLET** for improved malleability without using tiresome removable stylets
  - Improves ease of insertion
  - Increased value for money
- **COMPATIBLE** with all aqueous based contrast media.
- **AVAILABLE FROM OUR AUSTRALIA LOGISTICS PARTNER**

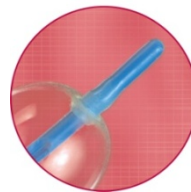

CLIFFORDHALLAMHEALTHCARE

1300 720 274 www.ch2.net.au

Ordering Information	Size	Code	CH2 Code
<b>Rocket HSG &amp; Sonohystogram Catheter</b> Length 280mm, balloon inflation syringe pre-attached. Maximum balloon volume: 3ml. Sterile, for single patient use, peel pouch packed in cartons of 10 catheters.	5FG	<b>R57422</b>	<b>2067145</b>
	7FG	<b>R57424</b>	<b>2067153</b>
<b>Rocket HSG &amp; Sonohystogram Catheter with integral wire stylet.</b> Length 280mm, balloon inflation syringe pre-attached. Maximum balloon volume: 3ml. Sterile, for single patient use, peel pouch packed in cartons of 10 catheters.	5FG	<b>R57422-00-ST</b>	<b>2067161</b>
	7FG	<b>R57424-00-ST</b>	<b>2067170</b>

