


HSG Procedure Packs

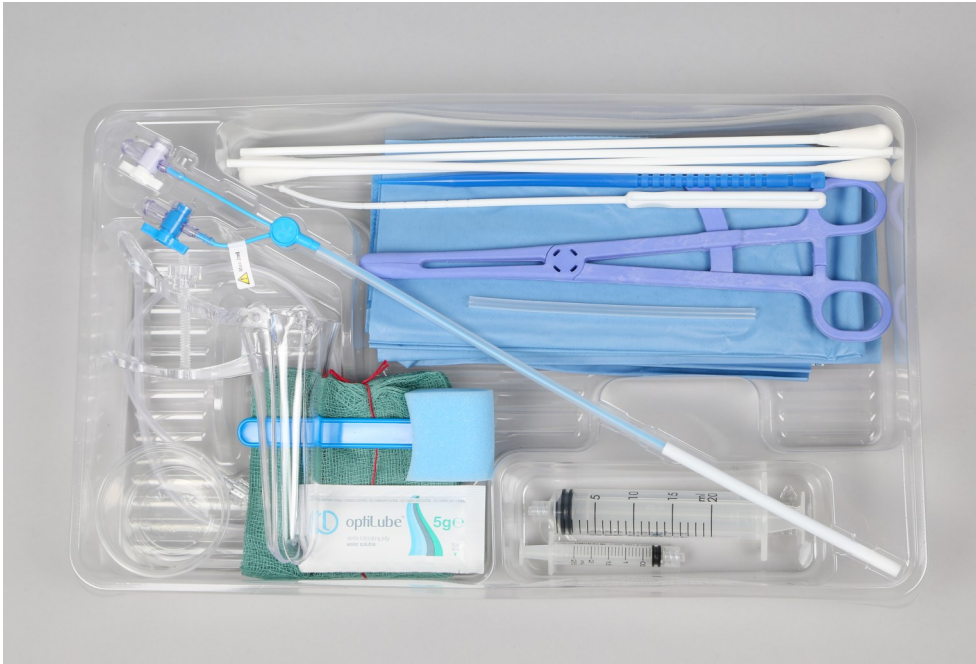

Rocket Medical HSG Procedure Packs provide all the necessary single use items for cost effective uterine canalisation for hysterosalpingography.

CONTENTS:

- HSG Catheter
- Vaginal speculum with screw retention
- Sponge stick for prep solution.
- Sponge holder with 5 x Green, X-ray detectable swabs
- Cervical dilator
- Water soluble lubricant
- 4 x 30cm swab sticks
- Filling kwill
- 20mL & 2.5mL syringes
- 30cm contrast extension tube

Rocket Medical HSG Catheters have these well-established and trusted features:

- **LATEX-FREE** balloon minimises risk of patient tissue reactions.
- **SMOOTH TIP AND LOW PROFILE BALLOON** ensures low trauma insertion and removal – ideal for the nulliparous cervix.
- **REMOVABLE STIFFENING SHEATH** to aid insertion in the heavily retroverted or antverted cervix.
- **CONTRAST PATH TAP** ensures retention of contrast when further volumes of media are required.
- **INTEGRAL WIRE STYLET** for improved malleability without using tiresome removable stylets (ST versions)
 - ◇ Improves ease of insertion
 - ◇ Increased value for money
- **COMPATIBLE** with aqueous based contrast media.

Description	Size	Order Code
HSG Procedure Pack HSG Catheter Length 280mm, Maximum balloon volume: 2ml. Sterile, for single use, in cartons of 10 procedure pack trays.	5FG	R57422-00-PK
	7FG	R57424-00-PK
HSG Procedure Pack with Stylet. HSG Catheter length 280mm, Maximum balloon volume: 2ml. Sterile, for single use, in cartons of 10 procedure pack trays	5FG	R57422-ST-PK
	7FG	R57424-ST-PK