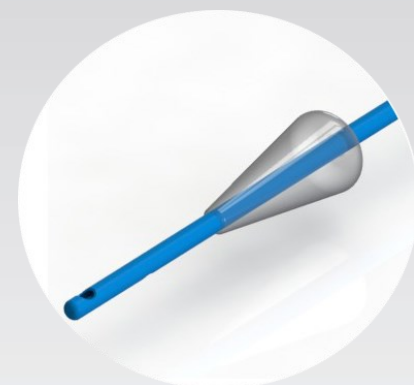

**The Rocket Medical Hysterocone™ Catheter has these features:**

- **MOVEABLE SILICONE CONE** ensures a soft, trauma free seal at the external os.
- **1cm GRADUATIONS** on the distal 7cm to allow the cone to be accurately positioned to match each patient's needs.
- **SMOOTH TIP AND 6FG TUBE** ensures simple, low trauma insertion, ideal for the nulliparous cervix.
- **REMOVABLE HOLLOW STIFFENING TUBE** to aid insertion, and to maintain rigidity during contrast injection.
- **CONTRAST PATH TAP** ensures retention of contrast when further volumes of media are required.
- **COMPATIBLE** with aqueous and oil –based x-ray and ultrasound contrast.

Description	Size	Order Code
<b>Hysterocone™ - 6Fg, HSG and HyCoSy Catheter</b> Active length 280mm, with 13mm cone, removable hollow stiffener tube and contrast path tap. Sterile, for single use, peel pouch and tray packed in cartons of 10.	6Fg	<b>R57430-CONE</b> <b>CH2: 2376491</b>

The Rocket Medical **Hysterocone™ Catheter** combines the features of both traditional cone and HSG catheters into a simple to use device for performing uterine canalisation for hysterosalpingography and hysterosalpingo contrast sonography.



CLIFFORDHALLAMHEALTHCARE

1300 720 274 www.ch2.net.au