


The Rocket Medical **Hysterocone™ Catheter** combines the features of both traditional cone and HSG catheters into a simple to use device for performing uterine canalisation for hysterosalpingography and hysterosalpingo contrast sonography.

The Rocket Medical Hysterocone™ Catheter has these features:

- **MOVEABLE SILICONE CONE** ensures a soft, trauma free seal at the external os.
- **1cm GRADUATIONS** on the distal 7cm to allow the cone to be accurately positioned to match each patient's needs.
- **SMOOTH TIP AND 6FG TUBE** ensures simple, low trauma insertion, ideal for the nulliparous cervix.
- **REMOVABLE HOLLOW STIFFENING TUBE** to aid insertion, and to maintain rigidity during contrast injection.
- **CONTRAST PATH TAP** ensures retention of contrast when further volumes of media are required.
- **COMPATIBLE** with aqueous and oil –based x-ray and ultrasound contrast media.

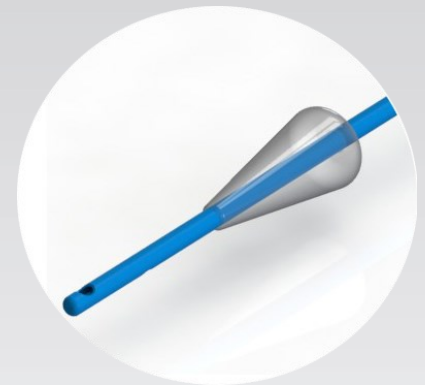

Description	Size	Order Code
Hysterocone™ - 6Fg HSG and HyCoSy Catheter Active length 280mm, with 13mm cone, removable hollow stiffener tube and contrast path tap. Sterile, for single use, peel pouch and tray packed in cartons of 10.	6Fg	R57430-CONE