


# HSG and HyCoSy Catheters

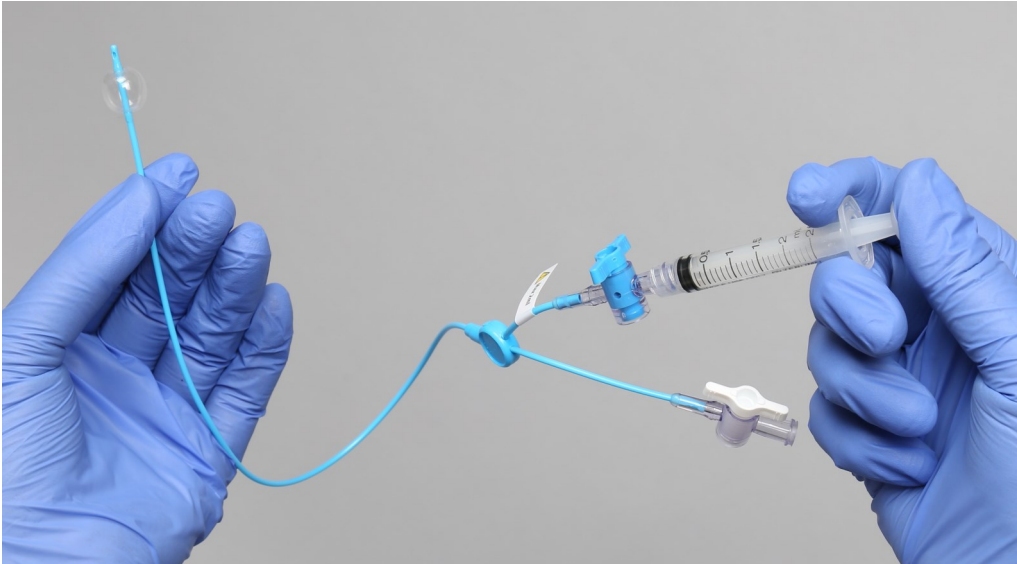

**Rocket Medical HSG Catheters provide safe and effective uterine canalisation for hysterosalpingography and hysterosalpingo contrast sonography.**

#### CONTENTS:

- HSG and HyCoSy Catheter
- Introducer sheath
- Pre-attached syringe
- Protective balloon cover

**Rocket Medical HSG and HyCoSy Catheters have these well-established and trusted features:**

- **LATEX-FREE** balloon anchors the catheter and provides a contrast seal.
- **SMOOTH TIP AND LOW PROFILE BALLOON** ensures low trauma insertion and removal – ideal for the nulliparous cervix.
- **REMOVABLE STIFFENING SHEATH** to aid insertion in the heavily retroverted or anteverted cervix.
- **CONTRAST PATH TAP** ensures retention of contrast when further volumes of media are required.
- **INTEGRAL WIRE STYLET** for improved malleability without using tiresome removable stylets (ST versions)
  - ◊ Improves ease of insertion
- **COMPATIBLE** with aqueous contrast media.

Description	Size	Order Code
<b>HSG and HyCoSy Catheter</b> HSG Catheter, length 280mm, Maximum balloon volume: 2ml. Sterile, for single use, in cartons of 10 units.	5FG	<b>R57422</b>
	7FG	<b>R57424</b>
<b>HSG and HyCoSy Catheter with Stylet</b> HSG Catheter, length 280mm. Maximum balloon volume: 2ml. Sterile, for single use, in cartons of 10 units.	5FG	<b>R57422-00-ST</b>
	7FG	<b>R57424-00-ST</b>